

# THE NEW DIMENSION IN PREVENTION

## PROFISIL® FLUORIDE VARNISH PLUS

FLUORIDE VARNISH

### A new dimension of varnish gel

**Profisil® Fluoride Varnish Plus** contains **5% sodium fluoride** in a mucosa-friendly, pleasantly flavored dimethicone gel. For optimal results, it contains soluble calcium and phosphate ions to promote uptake. This formula **is not made with ethyl alcohol or colophony (rosin)**, thus reducing the potential of soft tissue irritation or rosin-based allergies.

### Product Properties

- **Profisil® Fluoride Varnish Plus** contains soluble calcium and phosphate ions and can support the remineralization of tooth enamel.
- Due to its remineralising properties, hydroxyapatite is also important to tooth enamel and has good acid resistance. It is a main component of tooth enamel and is used in tooth pastes to repair, protect and strengthen the enamel.
- Developed on the basis of **Profisil®**, **Profisil® Plus** combines all the positive properties of the already well-known fluoride varnish **Profisil®**.


Calcium  
+ Phosphate  
Hydroxyapatite



### BENEFITS FOR THE DENTAL PROFESSIONAL



- **Profisil® Fluoride Varnish Plus** can be applied in any direction to all affected tooth surfaces, migrating even to difficult to access areas
- A prophylactic cleaning is not required prior to use
- Unnecessary to dry teeth prior to application of the varnish gel
- Lubricious, non-clumping formula  
» No need to replace gummy brushes
- Constant viscosity and easy to spread

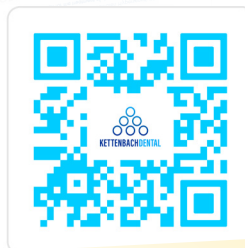

Hear Amanda Hill RDH on Profisil Plus!

# THE NEW DIMENSION IN PREVENTION PROFISIL® FLUORIDE VARNISH PLUS

FLUORIDE VARNISH



## BENEFITS FOR THE PATIENT

- Patients appreciate the smooth mouth feel and refreshing lubricating properties
- Protective coating to tooth surfaces by occluding tubules with calcium fluoride crystals
- No bitter taste or burning » not made with ethyl alcohol
- Non-allergenic dimethicone formula
- Suitable for all ages and for children with primary teeth
- This formula is designed to securely adhere to the tooth surface for several hours while releasing fluoride ions up to 24 hours
- Transparent and colorless after applying for optimal aesthetics

Calcium  
+ Phosphate  
Hydroxyapatite

Step-by-step to a  
prophy-clean feel!



Scan me

Creates a silky smooth surface and an instant  
pain-free, refreshing feeling in the mouth!



## Product range

The varnish is available in following versions:

MINT

BERRY

UNFLAVORED

Available as single dose unit  
incl. a brush for easy application.

<b>Normal pack</b>	Mint		14822	
	Berry	50 x 0.5 mL each	14824	\$102.50
	Unflavored		14826	
<b>Combi pack</b>		10 x 0.5 mL Mint		
		10 x 0.5 mL Berry	14821	\$68.50
		10 x 0.5 mL Unflavored		
<b>Bulk pack</b>	Mint		14827	
	Berry	300 x 0.5 mL each	14828	\$331.00
	Unflavored		14829	
<b>Superbulk pack</b>	Mint		14830	
	Berry	1000 x 0.5 mL each	14831	\$921.00
	Unflavored		14832	

Composition: Polydimethylsiloxane, 5% sodium fluoride (equivalent to 22.6 mg/g fluoride ions), Hydroxyapatite, soluble Phosphate salts, soluble Calcium salts, Steviol glycoside\*, xylitol\*, film-forming organic polymer, flavor\*.  
This formulation is free of gluten, peanuts, nuts, soy, sesame, corn, shellfish, eggs, or dairy products.

\* Not included in the unflavored version.



CALL 877-532-2123

Webshop: [www.kettenbachusa.com](http://www.kettenbachusa.com)

Website: [www.kettenbach-dental.us](http://www.kettenbach-dental.us)