




New Cartridge System



# SO ADHERENT IS **MUCOPREN®** SOFT

**A-SILICONE FOR SOFT RELININGS**

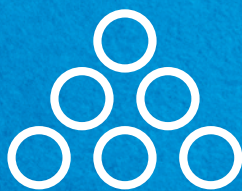

**KETTENBACHDENTAL**  
Simply intelligent



# MUCOPREN® SOFT

## Prosthetic relining

**Mucopren® Soft** is a soft long-term relining material based on A silicone for **direct and indirect application**.

### Special long-term adhesion:

Excellent adhesion to acrylic dentures through patented Mucopren® Adhesive

### Unique sealing through its very smooth, hydrophobic silicone surface:

Inhibits bacterial colonization thus ensuring hygienic conditions on the prosthesis<sup>1)</sup>



Application video

### Simple treatment & fast processing

Faster setting times and shorter waiting times ensure a fast and smooth process in the practice and laboratory.

### Supports the healing process of implants

Also suitable for use in implantology for trouble-free osseointegration.



### High wearing comfort

The definitive resilient reline is completely odorless and tasteless and makes wearing the prosthesis as pleasant as possible.



<sup>1)</sup> Result from in-vitro studies in which it was shown that the growth of *Candida albicans* was prevented thanks to the smooth silicone surface

# HOLDS TO ITS PROMISES



## FURTHER INDICATIONS



- For the creation of **permanently soft relinings**
- For relieving pressure on the surgical site after insertion of **implants**
- For balancing mucous membrane areas with different resilience
- For relieving sharp bony ridges
- For relining **interim prostheses**
- As **mucous membrane conditioner** in case of prosthetic stomatitis
- For creating and relining **obturators**

## Benefits of using the new adhesive

» ONLY A **SINGLE APPLICATION** NECESSARY



### CHAIRSIDE

3 MIN

The processing time is reduced to only **3 minutes** following mouth removal. No water bath necessary.



### LABSIDE

3 MIN

The processing time is reduced to only **3 minutes** in the water bath or pressure pot.

<sup>2)</sup> Refers to the optimized product version of Mucopren® Adhesive, first introduced in 2023.



# MUCOPREN® SOFT – THE RANGE

## Mucopren® Soft

<b>Mucopren® Soft Basic Set</b>	REF 28105	50 mL Mucopren® Soft, 50 mL Mucopren® silicone sealant, 10 mL Mucopren® Adhesive, 6 mixing tips white, 16 mixing tips orange, 1 brush holder, 20 single-use brushes, 1 steel bur, accessories
<b>Mucopren® Soft Normal pack</b>	REF 15687	2 x 50 mL Mucopren® Soft, 6 mixing tips white
<b>Mucopren® silicone sealant</b>	REF 15686	50 mL Mucopren® silicone sealant, 8 mixing tips orange
<b>Mucopren® Adhesive</b>	REF 14203	10 mL Mucopren® Adhesive



## Technical data

Product name	Mixing ratio	Working time at 73 °F ≤	Intraoral setting time at 95 °F ≥	Water bath	End of setting ≥	Shore hardness (approx.)
Mucopren® Soft – direct relining	1:1	2:15 min.	3:15 min.	Not applicable	5:30 min.*	A 28
Mucopren® Soft – indirect relining	1:1	2:15 min.	Not applicable	3:00 min.	> 3 min.	A 28
Mucopren® Soft silicone sealant	1:1	1:30 min.	Not applicable	Not applicable	6:30 min. at 23°C	A 33

\*Total setting time (removal time from mouth) from start of mixing

Kettenbach LP  
16052 Beach Blvd. Ste. 221  
Huntington Beach, CA 92647, USA  
Phone: 877-532-2123  
Fax: 888-892-9820  
Website: [www.kettenbach-dental.us](http://www.kettenbach-dental.us)  
Web shop: [www.kettenbachusa.com](http://www.kettenbachusa.com)  
Custservice@kettenbachusa.com

0210976\_2525

  
**KETTENBACHDENTAL**  
Simply intelligent